

BLICKMAN



Includes:

- Bassinet Kits
- Bassinet Frames
- Basket & Mattress

Bassinets

Transport baby safely and comfortably

OUR MISSION

Customer Focus

Exceed customer expectations through a commitment to continuous improvement of our products and services.

Quality Products

Offer the marketplace a broad range of high quality, cost effective, and durable products of superior craftsmanship.

Organizational Excellence

Strive for excellence in all we do by fostering a participative, team based environment, which nurtures the active involvement of all members of the organization.

Marketplace Recognition

Earn the respect, admiration, and loyalty of our people, customers, and suppliers by treating them with dignity and respect, embracing diversity, and developing and maintaining long term partnerships.

Collaborative Workplace


Encourage all employees to work safely and efficiently, contribute innovative ideas, and personally benefit from the ongoing success of the company.

The logo for Blickman, featuring the word "BLICKMAN" in a bold, white, sans-serif font. The text is contained within a dark blue rectangular box with a thin white border. The box has a slight 3D effect with a lighter blue shadow on the right side.

THE BLICKMAN DIFFERENCE

- Our bassinets are renowned for their reliable 304 stainless steel, welded construction, durability and optimal performance which has made them the preferred choice among maternity nurses for decades.
- Our inverted “U-frame” bassinet design provides a stable and comfortable resting place for newborns. Our standard 4” ball bearing swivel casters, with two brakes, ensure safe transport.
- We offer various storage choices such as shelves, drawers and enclosed cabinets, allowing you to customize the bassinet to suit your individual needs.
- Our clear acrylic infant basket includes a cardholder for baby ID and rests in specially formed support rods so it can be set in the flat or Trendelenburg position.

Key Product Features

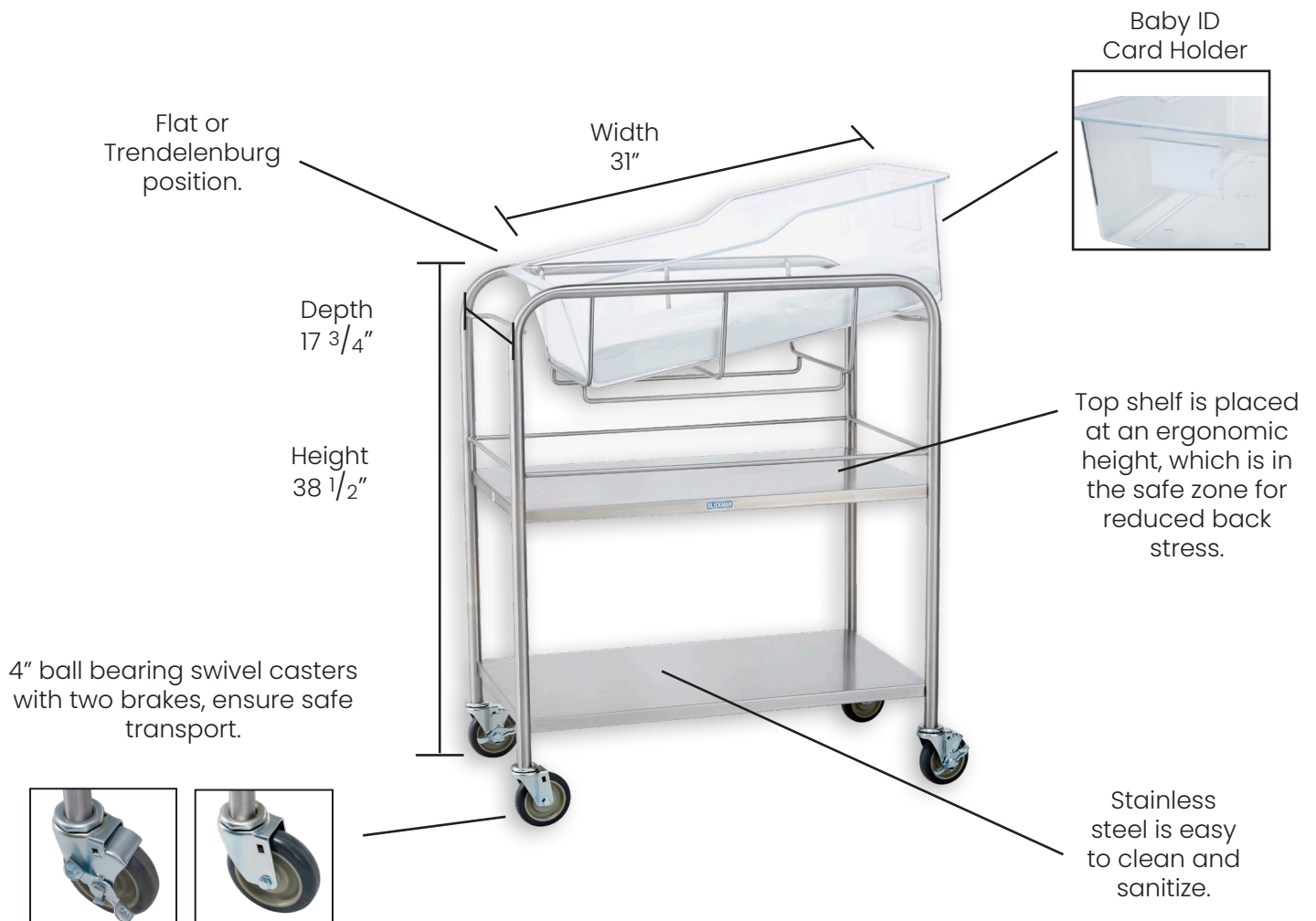

Bassinets

Bassinet Kits

We offer bassinet kits so you can get everything you need in one order. Our bassinet kits include:

- ✓ The stainless steel bassinet of your choosing
- ✓ The acrylic basket
- ✓ The foam mattress pad



**Now
Included**

BASSINET KITS

All of our bassinets are the same dimensions:

Width: 31" x Height: 41 1/2" x Depth: 17 3/4"



Description:

- Upper shelf has four sided guardrail to protect supplies.
- Two shelves are 8" and 23" from the floor.
- Shelf dimensions (WxHxD): 31" x 17 3/4"

Description:

- Storage drawer. Counter weighted for stability.
- One shelf, 8" from the floor.
- Drawer inside dimensions (WxHxD): 19 3/4" x 6 1/2" x 13 7/8"

Description:

- Two 17" deep storage compartments with doors.
- Right side has removable shelf.
- Cabinet inside dimensions (WxHxD): 29" x 13 3/4" x 15 7/8"

Description:

- Two open 17" deep storage compartments.
- Right side includes non-removable shelf.
- Cabinet inside dimensions (WxHxD): 29" x 13 3/4" x 16 5/8"

Name	Bassinet Kit with Open Shelves	Bassinet Kit with Drawer	Bassinet Kit with Closed Cabinet	Bassinet Kit with Open Cabinet
Model #	8044SS-KIT	8047SS-KIT	8048SS-KIT	8049SS-KIT
Item #	1118044KIT	1118047KIT	1118048KIT	1118049KIT
Shipping Weight	69 lbs.	98 lbs.	98 lbs.	77 lbs.

Custom Designs Available Upon Request

Please email us at info@blickman.com

Bassinets

Parts Ordered Separately

BASSINETS

All of our Bassinets are the same dimensions:

Width: 31" x Height: 38 1/2" x Depth: 17 3/4"



Name	Bassinet with Open Shelves	Bassinet with Drawer	Bassinet with Closed Cabinet	Bassinet with Open Cabinet
Model #	8044SS	8047SS	8048SS	8049SS
Item #	1118044000	1118047000	1118048000	1118049000
Ship Weight	53 lbs.	83 lbs.	82 lbs.	61 lbs.

ACCESSORIES



Basket

Dimensions: 31 1/8" x 9 7/8" x 16 1/8"

Seamless, clear acrylic bassinet basket with cardholder for easy identification.



Mattress

Dimensions: 26" x 2" x 12 1/2"

2" thick foam padding with stap check cover. Flame retardant.

Name	Basket	Mattress
Model #	3070P	3069MP
Item #	9841307001	8441306901
Ship Weight	12 lbs.	3 lbs.

Custom Designs Available Upon Request

Please email us at info@blickman.com



Bassinets

Construction Specifications

Construction Specifications

Universal Dimensions	Bassinets: Width: 31" x Height: 38 1/2" x Depth: 17 3/4" Bassinet Kits: Width: 31" x Height: 41 1/2" x Depth: 17 3/4"
Materials	304 stainless steel
Gauge	Legs: 16GA Shelves: 18 GA Body Parts and Removable Inside Shelves: 20GA
Welding	All welding by heliarc method and exposed welds are ground smooth.
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish.
Casters	4" ball bearing, swivel casters, two with brakes.
Notes	All models utilize specially formed support rods so basket can be set in the flat or Trendelenburg position. Infant basket and mattress pad can be sold separately.

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.

Custom Designs Available Upon Request
Please email us at info@blickman.com

BLICKMAN



OUR PROMISE

For over 100 years, Blickman has been known throughout the healthcare industry for exceptional product quality, and we pride ourselves on providing top notch and personalized service to ensure our customers are 100% satisfied. If you are ever not completely satisfied for any reason, we will make it right. And with our equipment backed by our Limited Lifetime Warranty, you can be certain you have made the very best choice for your healthcare facility.