

# BLICKMAN



**Includes:**

- Round Hampers
- Square Hampers
- Space Saver Hamper
- Hamper Bags

***Hampers***  
*For elimination of dirty linens*

# OUR MISSION

## **Customer Focus**

Exceed customer expectations through a commitment to continuous improvement of our products and services.

## **Quality Products**

Offer the marketplace a broad range of high quality, cost effective, and durable products of superior craftsmanship.

## **Organizational Excellence**

Strive for excellence in all we do by fostering a participative, team based environment, which nurtures the active involvement of all members of the organization.

## **Marketplace Recognition**

Earn the respect, admiration, and loyalty of our people, customers, and suppliers by treating them with dignity and respect, embracing diversity, and developing and maintaining long term partnerships.

## **Collaborative Workplace**

Encourage all employees to work safely and efficiently, contribute innovative ideas, and personally benefit from the ongoing success of the company.


The logo for Blickman, featuring the word "BLICKMAN" in a bold, white, sans-serif font. The text is contained within a dark blue rectangular box with a thin white border. The box has a slight 3D effect with a lighter blue shadow on the right side.

# THE BLICKMAN DIFFERENCE

- Our hampers facilitate easy and safe disposal of contaminated linens. They all move easily over floors on 2" or 3" casters. Our hampers have proven to be stable when empty or full due to their wide bases and added counterweights on covered models. A bottom brace holds and protects the hamper bag from contact with the floor.
- Our hampers come in either chrome or stainless steel and we offer different shapes and sizes to fit your needs and floor plans. Additionally, we offer covered models to control airborne pathogens. All covered hampers are foot operated to help prevent cross contamination.
- For MR conditional environments, we offer two MR conditional models: one covered and one uncovered, for your choosing.
- We offer two optional hamper bags, nylon and water soluble, that fit our round hampers.

## *Key Product Features*

---

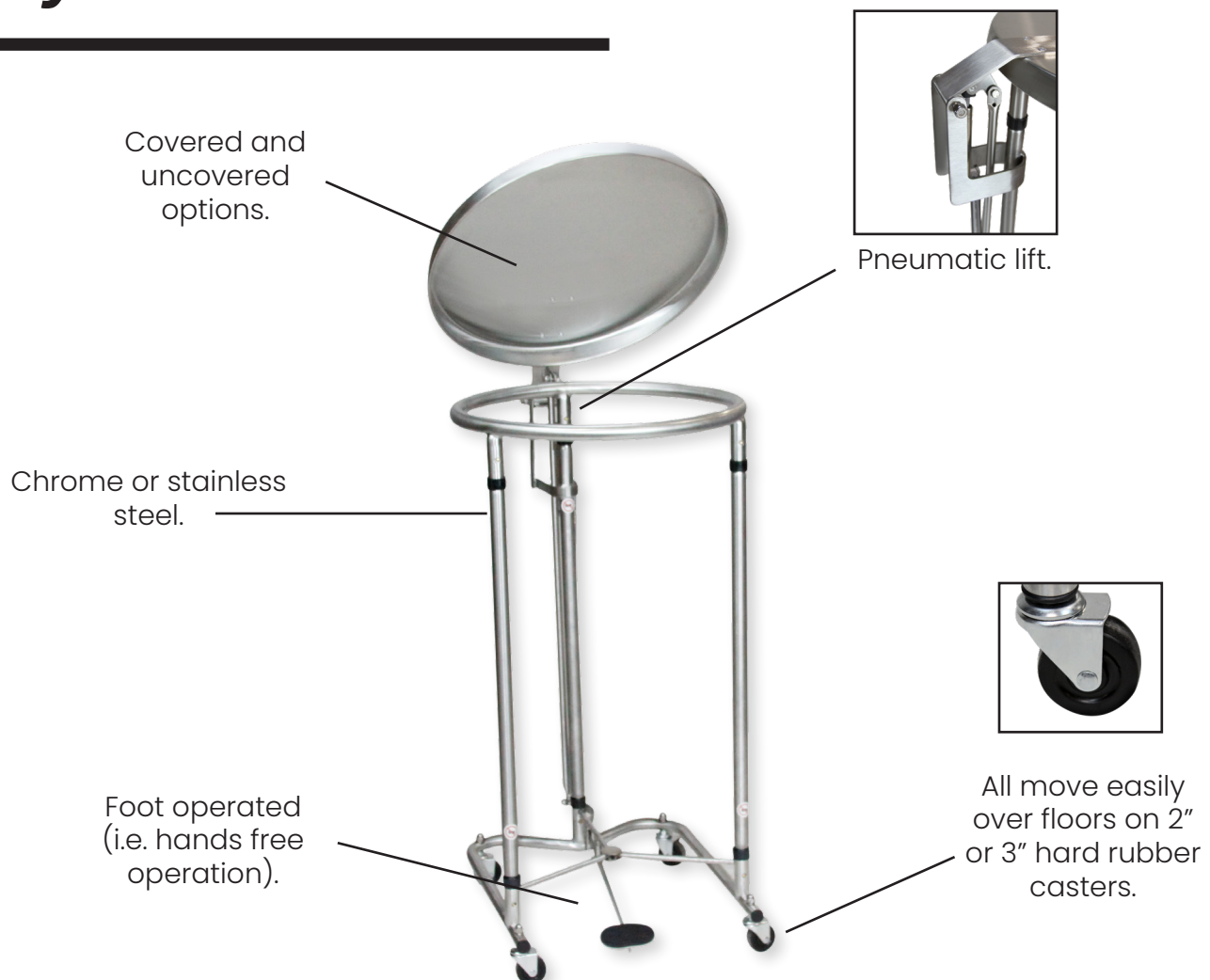


# Hampers

Round

## ROUND HAMPERS

**Description:** We offer covered and uncovered round hampers in both stainless steel and chrome. They are 35" tall and roll on 2" hard rubber swivel casters. They come in two sizes 18" diameter and 25" diameter. The covered models operate via a spring-loaded lift mechanism that opens and closes the stainless steel lid. The pneumatic top allows for pin drop quiet operation.

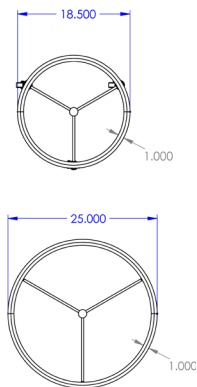

Name	Uncovered		Covered	
Lid	Chrome	Stainless Steel	Chrome	Stainless Steel
<b>Bag Size</b>	<b>18"</b>			
Model #	2018	7774SS	2118	7774SS-LF
Item #	0962018000	0927774000	096211800P	092777420P
Shipping Weight	11 lbs.	11 lbs.	21 lbs.	25 lbs.
<b>Bag Size</b>	<b>25"</b>			
Model #	2025	8773SS	N/A	8773SS-LF
Item #	0962025000	0928773000	N/A	092877320P
Shipping Weight	12 lbs.	13 lbs.	N/A	35 lbs.

**Hamper Bags:** 18" and 25" white nylon hamper bags and water soluble hamper bags

MR options available, see MR section for more details.



# Hampers

Square

## FOLDING ADJUSTABLE HAMPER

**Description:** Our folding chrome hamper adjusts to accommodate three hamper bag sizes. A bottom frame supports hamper bag from underneath. Travels on 2" hard rubber swivel casters.



**Dimensions:**

- **18" bag size:** 19" x 36 1/2" x 18"
- **22" bag size:** 19" x 40 1/2" x 22"
- **25" bag size:** 19" x 42 1/4" x 25"

Name	Folding Adjustable Hamper
Model #	2230
Item #	0962230000
Shipping Weight	18 lbs.

## SECURE CARE HAMPER

**Description:** Our secure care chrome hamper supports 36 - 42 gallon multipurpose bags. Its' polycoated steel lid opens using foot operated lift mechanism. Includes hamper storage pouch. Travels on 3" hard rubber swivels casters.



**Dimensions:**

- 18 5/8" x 37 3/4" x 19 1/2"

Name	Secure Care Hamper
Model #	2010
Item #	0962010000
Shipping Weight	26 lbs.



# Hampers

Triangular

## SPACE SAVER HAMPER

**Description:** Our space saver hamper offers hands free loading and unloading with a convenient non-slip foot pedal. Sturdy all stainless steel construction with bumper guards on lid and legs to protect equipment and walls. The solid bottom shelf provides extra security for heavy loads and a spring-loaded lift mechanism provides a pin drop pneumatic top for quiet operation. Travels on 3" hard rubber swivel casters.



Name	Space Saver Hamper
Model #	230ISS
Item #	096230100P
Shipping Weight	24 lbs.

# Hampers

Accessories

## HAMPER BAGS



**Disposable Water Soluble  
Plastic Hamper Bag**



**Reusable  
Nylon Hamper Bag**



**Secure Care  
Hamper Pouch**

Part #	Description	Shipping Weight
98390000WS	100 Plastic Water Soluble Hamper Bags	9 lbs.
9839777201	Bag, 18" Diameter Ring Hamper Nylon White	1 lb.
9839777302	Bag, 25" Diameter Ring Hamper Nylon White	1 lb.
0962012000	Secure Care Hamper Pouch for Labels	1 lb.





# Hampers

Construction  
Specifications

## Construction Specifications

Universal Dimensions	<b>Hamper Bag Sizes:</b> 18", 22", 25", 26"
Materials	304 stainless steel Chrome
Gauge	<b>Stainless Steel:</b> 1" diameter, 18GA tubing <b>Chrome:</b> 1" diameter tubing
Welding	All welding by heliarc method and exposed welds are ground smooth.
Finish	Exposed surfaces are free from pit marks, weld seams, scale and are polished with #4 polish to a satin finish.
Casters	2" & 3" hard rubber swivel casters (some with locks).
Notes	MR conditional hampers available where indicated.

## Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at [info@blickman.com](mailto:info@blickman.com).



**BLICKMAN**



## OUR PROMISE

For over 100 years, Blickman has been known throughout the healthcare industry for exceptional product quality, and we pride ourselves on providing top notch and personalized service to ensure our customers are 100% satisfied. If you are ever not completely satisfied for any reason, we will make it right. And with our equipment backed by our Limited Lifetime Warranty, you can be certain you have made the very best choice for your healthcare facility.