

BLICKMAN



Utility Carts

Includes:

- Light Duty
- Medium Duty
- Heavy Duty

Convenient transport for medical equipment

OUR MISSION

Customer Focus

Exceed customer expectations through a commitment to continuous improvement of our products and services.

Quality Products

Offer the marketplace a broad range of high quality, cost effective, and durable products of superior craftsmanship.

Organizational Excellence

Strive for excellence in all we do by fostering a participative, team based environment, which nurtures the active involvement of all members of the organization.

Marketplace Recognition

Earn the respect, admiration, and loyalty of our people, customers, and suppliers by treating them with dignity and respect, embracing diversity, and developing and maintaining long term partnerships.

Collaborative Workplace


Encourage all employees to work safely and efficiently, contribute innovative ideas, and personally benefit from the ongoing success of the company.

The logo for Blickman, featuring the word "BLICKMAN" in a bold, white, sans-serif font. The text is contained within a dark blue rectangular box with a thin white border. The box has a slight 3D effect with a lighter blue shadow on the right side.

THE BLICKMAN DIFFERENCE

- Our utility carts are ideal for transporting medical equipment which need to be used in many different areas in the facility. We offer light duty, medium duty, and heavy duty options to suit all your needs.
- Common features include a push handle, open construction for easy access, guardrails for securing contents, and 4" swivel casters allow for easy maneuverability.
- They contain two or three solid stainless steel shelves providing ample storage space.
- The stainless steel construction allows for easy cleaning and sanitizing, while resisting corrosion and stains.
- We also offer a folding utility cart with wire baskets for any application such as movement of mail, patient records, and laboratory materials.

Key Product Features

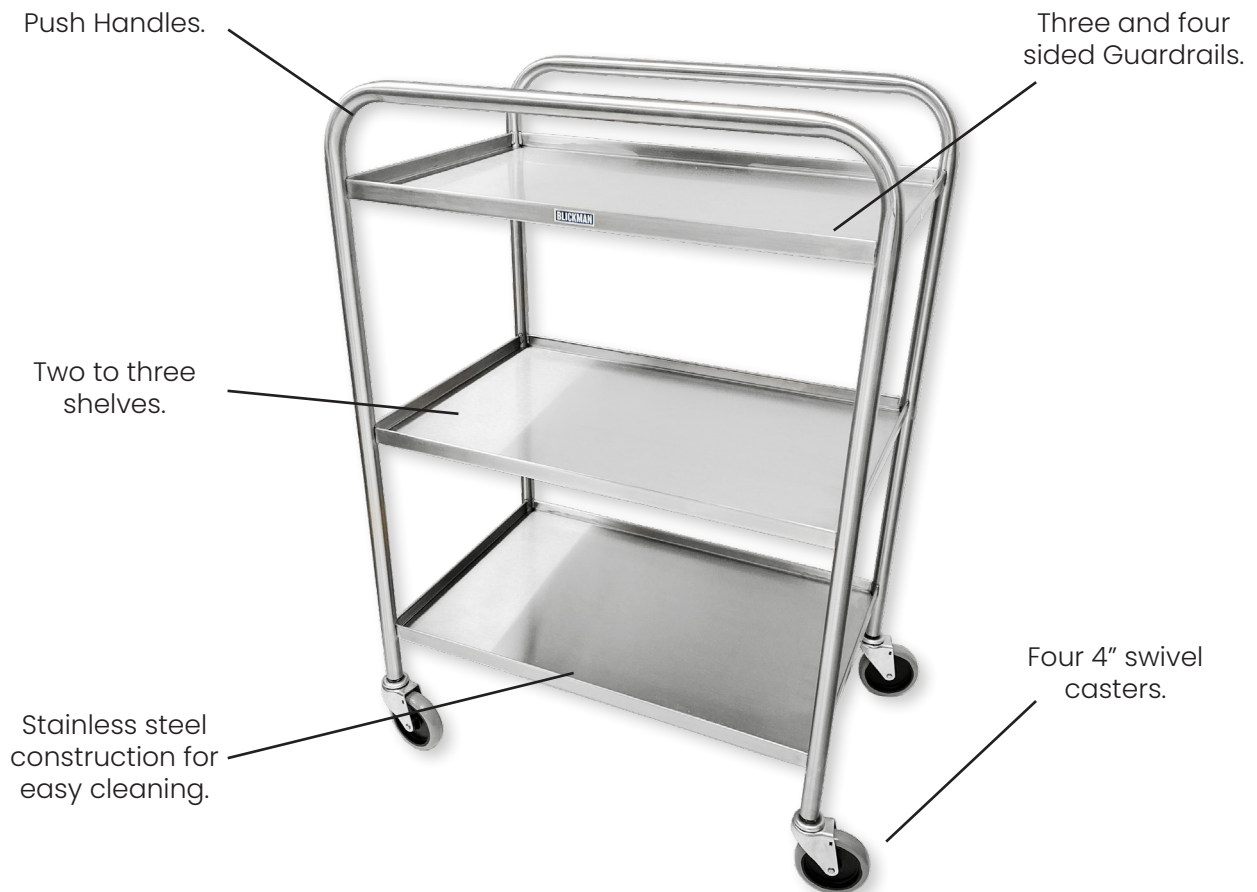

Utility Carts

Light Duty

FOLDING UTILITY CART

Description: Our folding utility cart is versatile enough for any application such as movement of mail, patient records, and laboratory materials. Choose any combination of 6" or 12" deep wire baskets to complete your cart. The metal frame is designed for easy assembly and folds for storage. Rubber stoppers support and protect the wire racks and provide a quiet ride on four, 3" hard rubber casters.

Folding Utility Cart Dimensions:

Width: 29 1/4" x Height: 38 1/2" x Depth: 19 1/2"



6" Wire Basket (2441)

Item #: 2422441000



12" Wire Basket (2442)

Item #: 2422442000

Name	Model #	Item #	Shipping Weight
Folding Utility Cart	2440	2422440000	21 lbs.
ORDER SEPARATELY			
6" Wire Basket	2441	2422441000	5 lbs.
12" Wire Basket	2442	2422442000	6 lbs.

Utility Carts

Light Duty

LIGHT DUTY CARTS

Description: Our light duty 7534 utility cart comes standard with three shelves and a push handle. The top two shelves have three sides and an open front to slide in equipment and the bottom shelf has no sides for easy access. We offer optional designs with additional attachments for differing needs.

Our light duty carts have the same dimensions:

Width: 28 1/4" x Height: 35 3/4" x Depth: 16"
Added IV pole: Width 30", Height: 67 1/2" - 109 1/2"



Top Shelf: Three Sided Shelf
Middle Shelf: Three Sided Shelf
Bottom Shelf: Open Shelf

Top Shelf: Three Sided Shelf
Middle Shelf: Three Sided Shelf
Bottom Shelf: Open Shelf

Top Shelf: Three Sided Shelf
Middle Shelf: Three Sided Shelf with Front Guardrail
Bottom Shelf: Shelf with Four Guardrails

Name	7534SS	7534 with IV pole	7534SS with Guardrails
Model #	7534SS	7534SS-IV	7534SS-GR
Item #	2427534002	2427534001	2427534000
Shipping Weight	40 lbs.	51 lbs.	42 lbs.

Custom Designs Available Upon Request
Please email us at info@blickman.com

Utility Carts

Medium/Heavy Duty

MEDIUM DUTY CARTS

Description: Our medium duty utility carts each come with three solid stainless steel shelves and travel on 4" swivel casters. Choose between our design that includes four sided shelves and width-wise handles or our design that includes open shelving and depth-wise handles.



Dimensions: 27 1/2" x 40" x 18 1/8"



Dimensions: 25" x 16 1/2" x 35 3/4"

Name	Side Rails Cart	End Rails Cart
Model #	7535SS	7536SS
Item #	2427535000	2427536000
Shipping Weight	65 lbs.	43 lbs.

HEAVY DUTY CART

Description: Our heavy duty utility cart comes with two solid stainless steel shelves with four sides each and travels on 5" swivel casters. Easy maneuverability with angled push handles.



Dimensions: 41 3/4" x 40" x 19 1/4"

Name	Heavy Duty Cart
Model #	7537SS
Item #	2427537000
Shipping Weight	59 lbs.

Custom Designs Available Upon Request
Please email us at info@blickman.com

Utility Carts

Accessories

IV POLE ATTACHMENTS

Description: Our IV pole attachment can fit any Blickman utility cart.



Name	7534SS with IV pole attachment	7535SS with IV pole attachment	7536SS with IV pole attachment	7537SS with IV pole attachment
Model #	7534SS-IV	7535SS-IV	7536SS-IV	7537SS-IV
Item #	2427534001	2427535001	2427536001	2427537001
Shipping Weight	51 lbs.	76 lbs.	54 lbs.	70 lbs.

GUARDRAILS & SIDING

Contact Blickman Customer Service about adding additional siding or guardrails.



Open Shelf



Lipped Shelf



Shelf with 4 Guardrails



3 Sided Shelf

Custom Designs Available Upon Request

Please email us at info@blickman.com

Utility Carts

Construction Specifications

Construction Specifications

Universal
Dimensions

Light Duty: Width: 28 1/4" x Height: 35 3/4" x Depth: 16"

Materials

Stainless steel
Rubber bumper handles on 7534SS Carts

Gauge

Light Duty: 20GA & 18GA
Medium/Heavy Duty: 16GA

Welding

All welding by spot and heliarc methods and exposed welds are ground smooth.

Finish

Exposed surfaces are free from pit marks, weld seams, scale and are polished with a #4 polish to a satin finish.

Casters

4" Swivel Casters

Shipping Information

All shipping weights are estimates and may be changed without prior notice. Specifications may be improved and dimensions changed without prior notice. Please contact customer service for estimated shipping costs at info@blickman.com.

Custom Designs Available Upon Request
Please email us at info@blickman.com

BLICKMAN



OUR PROMISE

For over 100 years, Blickman has been known throughout the healthcare industry for exceptional product quality, and we pride ourselves on providing top notch and personalized service to ensure our customers are 100% satisfied. If you are ever not completely satisfied for any reason, we will make it right. And with our equipment backed by our Limited Lifetime Warranty, you can be certain you have made the very best choice for your healthcare facility.