

Be Independent.™



SECURITY POLE & CURVE GRAB BAR™

BECAUSE YOUR GRAB BAR CAN GO ANYWHERE IN THE HOUSE




Point of Purchase
Also Available #1111

PRODUCT #1100



WATCH A DEMO



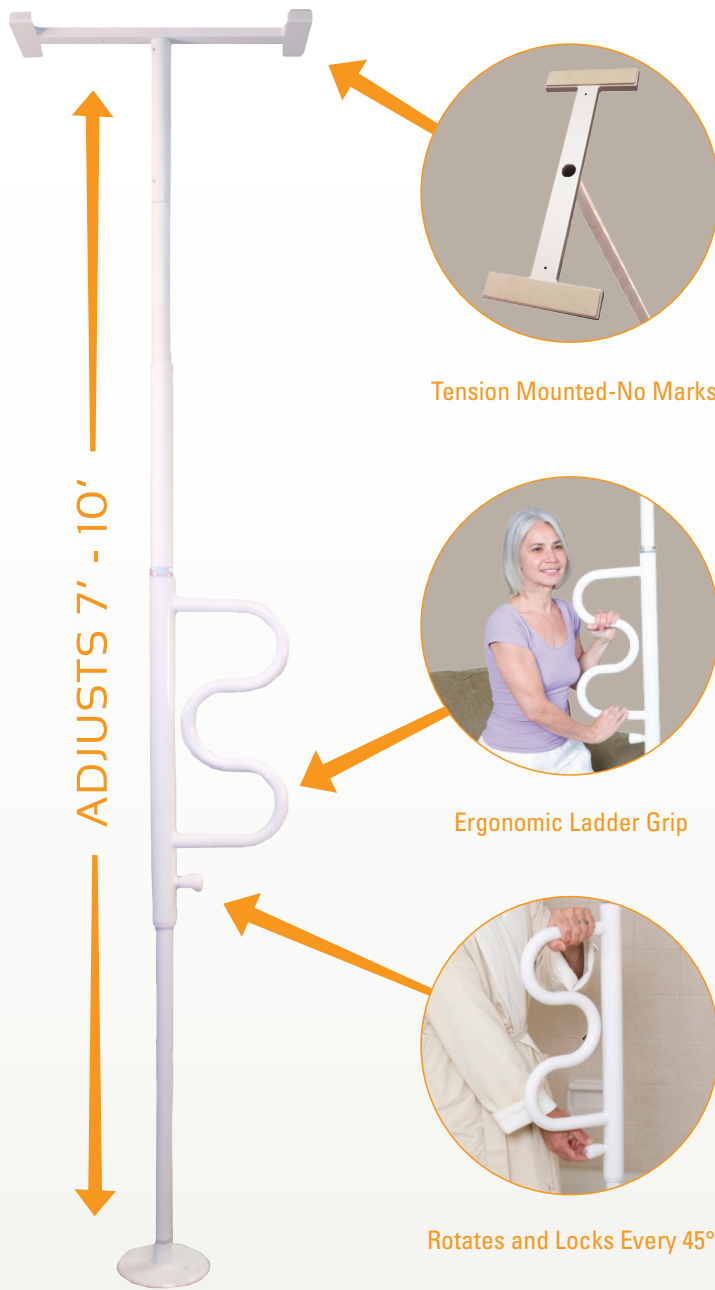

FEATURES & BENEFITS

- **Two Products in One** - Floor to ceiling pole with pivoting curve grab bar that locks in place every 45°.
- **Adjustable** - Fits ceiling heights 7'–10' ft.
- **Portable** - Main pole easily separates into two 53" pieces making it easy to transport.
- **Adaptable** - No need for wall mounts; can be installed in virtually any area of a room exactly where it is needed.

SPECIFICATIONS

Materials: Rust resistant zinc plated steel with durable powder coat finish

Fits Ceiling Heights: Range of 7'-10' ft

Weight Capacity: 300 lbs

Weight of Product: 20 lbs

Package Dimensions: 54.00" L x 12.50" W x 3.25" H; 22 lbs

Case of 2 Dimensions: 54.00" L x 12.50" W x 6.5" H; 45 lbs

COMES IN TWO COLORS

■ BLACK #1100-B

□ WHITE #1100-W

EASY ASSEMBLY (No Additional Tools Required)



Attach Top and Bottom Pieces Using Included Pins



Slide the Two 53" Pieces Together



Use Included Wrench to Tension Mount Between Floor and Ceiling