

Be Independent.™



# STABLE RAIL™

BECAUSE STABILITY MAKES STANDING EASY



PRODUCT #5800



WATCH A DEMO





Accessible Handle  
Provides Support



Includes a Bonus  
Organizer Pouch



Adjustable Legs  
Provide Stability



Accessorize the Stable Rail  
with the Stander Arm  
#5810



## FEATURES & BENEFITS

- **Swing Out Mobility Arm** – Pivots outward to provide improved walking stability.
- **Extendable Legs** – Provides added stability when standing up.
- **Ergonomic Cushion Handle** – Allows for easy transfer in and out of bed.
- **4-Pocket Organizer** – Provides storage space for handy items.
- **Anti-Slip Grips** – Secures rail in-between mattress and bed frame.
- **Height Adjustment** – Accommodates any home or hospital bed.
- **Easy Installation** – Installs in seconds with no tools required.

## SPECIFICATIONS

*Height Adjustment From Floor to Top of Handle: 21" - 44"*

*Height Adjustment From Mattress Base to Top of Handle: 17.5" - 22.5"*

*Rail Width: 19"*

*Height Range From the Floor to Horizontal Support Structure Under Mattress: 13.5" - 21.5"*

*Depth of Horizontal Support Structure Under Mattress: 25"*

*Weight Capacity: 300 lbs*

*Weight of Product: 10.5 lbs*

*Made of: Steel frame with durable powder coat finish*

*Package Dimensions: 28" L x 20.25" W x 2" H*

*Case of 4 Dimensions: 28" L x 20.25 W x 5.75" H; 32 lbs*


## EASY ASSEMBLY (No Additional Tools Required)




Adjust Base Tubes to  
Accommodate Left or  
Right Side of Bed



Attach Base Tubes  
to Main Handle

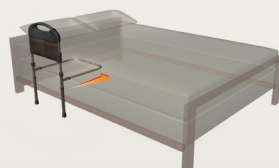

Insert Base Tubes Between  
Mattress and Box Spring

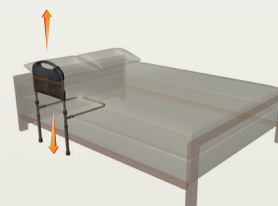

Adjust Height of Handle  
and Legs as Desired